

	FINISHED PRODUCT SPECIFICATION		Page No.	1 of 2
	Quality Control Department			
Product Name	NOVOMIX 110228			
Colour	Light Yellow			
Specification No.	FP/NOVOMPS/110228/01	Supersedes No.	NA	
Reference	In House	Effective Date	22/02/2024	
Ref. Annexure No.	A/SOP/QC/030/04	Review Month	Jan ' 2029.	

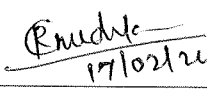
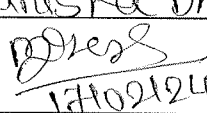
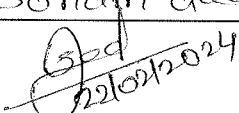
GENERAL INFORMATION

Pharmacopoeial reference : In House

Quantity to be sampled :	Analysis Sample	Control sample
	25g	75g

Shelf Life : 36 Months

Storage : Store in dry place, in well closed container, away from direct heat and light.

	PREPARED BY	REVIEWED BY	APPROVED BY
Name	Rajanikant Dudhkar	Manisha Dhere	Sonam Gao
Sign/Date	 17/02/24	 17/02/24	 22/02/2024
Designation	Sr. Q.C. Executive	Sr. Executive	Sr. Manager
Department	Quality control	Quality Control	Quality Assurance

	FINISHED PRODUCT SPECIFICATION		Page No.	2 of 2
	Quality Control Department			
Product Name	NOVOMIX 110228			
Specification No.	FP/NOVOMPS/110228/01			

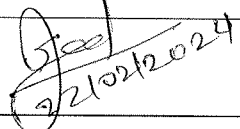
Qualitative Formula:

1. Hydroxypropyl methyl cellulose IP/USP
2. Polyethylene Glycol IP/USP
3. Talc IP/USP
4. Titanium Dioxide IP/USP
5. Red Iron Oxide USP
6. Yellow Iron Oxide USP

Sr. No.	Test	Specification	Reference
1.	Appearance	Light Yellow Powder.	GTP/QC/001
2.	Identification by IR	IR spectrum of test sample should be concordant with the IR spectrum of standard.	GTP/QC/009
3.	Colour Shade		GTP/QC/008
	1) Visually	Colour shade should match visually with Standard	
	2) DE (By colour spectrophotometer)	NMT 2.0	
4.	Particle Size (Done on wet slurry)	Not less than 99.00 % w/w passes through 100 mesh	GTP/QC/006
5.	Tapped Density	Not more than 1.2 g/ml	GTP/QC/003
6.	Ash Content	Not more than 30 % w/w	GTP/QC/004

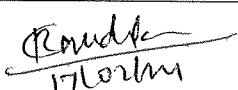
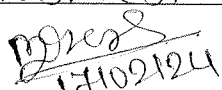
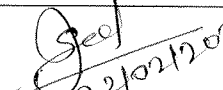
REVISION HISTORY FOR SPECIFICATION:

Revision No.	Revision Description	Effective Date
01	New Document	22/02/2024

	PREPARED BY	REVIEWED BY	APPROVED BY
Sign/Date	 17/02/24	 17/02/24	 22/02/2024

	FINISHED PRODUCT STANDARD TEST PROCEDURE (STP)		Page No.	1 of 4
	Quality Control Department			
Product Name	NOVOMIX 110228			
Colour	Light Yellow			
STP No.	FP/NOVOMSTP/ 110228/01			
Reference	In House	Supersedes No.	NA	
Ref. SPEC No.	FP/NOVOMPS/110228/01	Effective Date	22/02/2024	
Ref. Annexure No.	A/SOP/QC/030/08	Review Month	Jan 2024	

01.	Appearance : Instrument: Balance Take about 2 to 3 g of test sample and transferred into a clean and dry petri dish. Spread it uniformly by spatula and observe. Specification: Light Yellow Powder.
	02. Identification by IR: Instrument : Balance , IR Spectrophotometer, Oven Take about 1 mg of sample to be analyzed and 100 mg IR grade KBr onto piece of weighing paper. The KBr should be store in the desiccator when not in used. Transfer the sample /KBr mixture to the agate mortar. Thoroughly disperse the mixture using the spatula or spoon. The pestle is then used with a grinding motion to homogenously disperse the mixture component. Completed homogeneity takes approximately 1-2 minutes of grinding the components with the pestle. Prepare a KBr pellet using an appropriate apparatus (Press Kit, E- Z press etc.) carefully examine the pallet film to ensure that it is intact and somewhat transparent. Run the pallet on suitable spectrophotometer from 4000 cm^{-1} to 650 cm^{-1} according to the operating instruments. All peak should be match the reference with relation to wavelength location and approximate magnitude location. Overlap the sample spectra with reference spectra. Note: Sample and KBr to be dried at 105°C for about 30 minutes before use Or By ATR Accessory:

	PREPARED BY	REVIEWED BY	APPROVED BY
Name	Rajankant Soudhakar	Manisha Dhere	Sonam Geel
Sign/Date	 17/02/24	 17/02/24	 22/02/2024
Designation	Sr. Q.C. Executive	Sr. Executive	Sr. Manager
Department	Quality control	Quality Control	Quality Assurance

2 of 4

Quality Control Department

NOVOMIX 110228

FP/NOVOMSTP/ 110228/01

Specification: IR spectrum of test sample should be concordant with the IR spectrum of standard.

Note: Dispersion media and sample quantity may vary as per requirement but final concentration should remain same.

Sign/Date

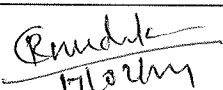
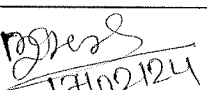
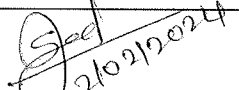
Boudk
17/02/24

17102124

~~Sec~~
~~2402/2024~~

	FINISHED PRODUCT STANDARD TEST PROCEDURE (STP)		Page No.	3 of 4
	Quality Control Department			
Product Name	NOVOMIX 110228			
STP No.	FP/NOVOMSTP/ 110228/01			

	Compare the colour shade cards with a standard reference colour shade card. Specification: Colour shade should match visually with Standard. 2) DE (By colour spectrophotometer): Compare the colour of sample shade cards with a reference standard already stored in the colour matching system. Specification: NMT 2.0
04.	Particle Size (Done on wet slurry): Instrument: Balance, Stirrer, Hot air oven Prepare about 5% dispersion (5 g in 95 g $\pm$ 1 g of distilled water) and stir the dispersion for about 45 minutes with constant stirring. Pour this dispersion on the top of the 100 mesh tarred mesh in a form of fine stream. Sieve the remaining particles on the mesh without applying any pressure with the help of running water. Continue the wet sieving process until the emerging liquid appears free of particles. Remove the mesh and dry it at 105°C. Allow the mesh to come to room temperature in a desiccator before weighing. Repeat the operation until two successive weighing do not differ by more than 1 mg. Determine the weight of dried material on the mesh. Calculation: Percentage w/w passed through 100 mesh = $\frac{100 \times W1 - W5}{W1}$ Where, W1 = Weight of sample in g taken for mesh test. W2 = Weight of empty mesh W3 = Weight of mesh plus residue after drying W4 = Weight of mesh plus residue after drying (constant weight) W5 = Weight of material remained on the mesh after drying to constant weight = W4 - W2 Specification: Not less than 99.00 % w/w passes through 100 mesh.

	PREPARED BY	REVIEWED BY	APPROVED BY
Sign/Date		 17/02/24	 22/02/2024

	FINISHED PRODUCT STANDARD TEST PROCEDURE (STP)		Page No.	4 of 4
	Quality Control Department			
Product Name	NOVOMIX 110228			
STP No.	FP/NOVOMSTP/ 110228/01			

05. Tapped Density:

Instrument: Balance, Densitometer

Weigh 10 g of sample. Pour it in a 50 ml Graduated measuring cylinder. Set 100 taps for the densitometer. Note the volume (A) occupied by 10 g of the powder after 100 taps on the graduated cylinder.

Calculations: Tapped Density in g/ml = $\frac{\text{Weight of sample}}{\text{Volume A}}$

Specification : Not more than 1.2 g/ml

06. Ash Content :

Instrument: Balance, Muffle Furnace

Heat a silica crucible at 600 ° C ± 25° C for 15 minutes, allow to cool in a desiccator and weigh.

Take 1-2 g of sample to the crucible and weigh the crucible. Ignite the sample at 600 ° C ± 25° C for 45 minutes. Allow the silica crucible to come to room temperature in a desiccator before weighing.

Repeat the operation until two successive weighing do not differ by more than 0.5 mg

Calculation:

$$\% \text{ Ash} = \frac{(W4-W1)}{(W2-W1)} \times 100$$

Where, W1 = Weight of the empty crucible (g)

W2 = Weight of the crucible with test sample before ignition (g)

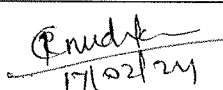
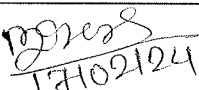
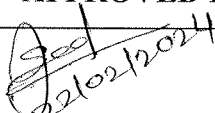
W3 = Weight of the crucible with residue after ignition (g)

W4= Weight of the crucible with residue after ignition (g) (Constant weight)

Specification: Not more than 30 % w/w.

REVISION HISTORY FOR STANDARD TEST PROCEDURE:

Revision No.	Revision Description	Effective Date
01	New Document	22/02/2024

	PREPARED BY	REVIEWED BY	APPROVED BY
Sign/Date	 17/02/24	 17/02/24	 22/02/2024